

TECHNICAL SHEET

Cage and Tray Specifications and Dimensions

Animal Model	Cat. N.	Type	Sizes (W x D x H) mm/inch	Area (W x D x H) mm/inch	Weight Kg
Multifloor Pen System	 4421X009	*X* type cage (compatible with tray Cat. N. 4401S389)	716 x 716 x 114 28.19 x 28.19 x 4.48	653 x 653 x 95 25.71 x 25.71 x 3.74	2.8
	 4401S389	Tray (compatible with *X* type cage Cat. N. 4421X009)	716 x 716 x 46 28.19 x 28.19 x 1.81	656 x 656 x 38 25.83 x 25.83 x 1.496	2.4
	 4541P009	*P* type cage (compatible with tray Cat. N. 4501K389)	870 x 740 x 130 34.25 x 29.13 x 5.118	797 x 679 x 112 31.38 x 26.73 x 4.409	3.8
	 4501K389	Tray (compatible with *P* type cage Cat. N. 4541P009)	866 x 786 x 46 34.09 x 30.94 x 1.81	800 x 720 x 40 31.5 x 28.35 x 1.575	3.6
Ferret Cage	 4541P008	*P* type cage (compatible with tray Cat. N. 4501K389)	870 x 740 x 130 34.25 x 29.13 x 5.118	797 x 679 x 122 31.38 x 26.73 x 4.803	3.8
Plastic Guinea Pigs Cage	 3020W008	*W* type solid base cage (compatible with tray Cat. N. 3020W389)	814 x 610 x 260 32.05 x 24.02 x 10.24	743 x 543 x 250 29.25 x 21.38 x 9.843	4.3
	 3020W009	*W* type perforated base cage (compatible with tray Cat. N. 3020W389)	814 x 610 x 260 32.05 x 24.02 x 10.24	743 x 543 x 250 29.25 x 21.38 x 9.843	4.3
	 3020W389	Tray (compatible with *W* type perforated base cage Cat. N. 3020W009)	809 x 609 x 46 31.85 x 23.98 x 1.81	749 x 549 x 40 29.49 x 21.61 x 1.575	2.4

Rack Specifications and Dimensions

Animal Model	Complete Rack Unit Cat. N.	Complete Rack unit Autovatering System Cat. N.	N. of Cages	Internal Height mm/inch	Overall Rack Dimensions (W x L x H) - mm/inch
Multifloor Pen System	 4X03B700CP	4X03B701CP	3	450 17.72	785 x 840 x 1840 30.91 x 33.07 x 72.44
	 4X06B700CP	4X06B701CP	6	450 17.72	1550 x 840 x 1840 61.02 x 33.07 x 72.44
	 4P02B700CP	4P02B701CP	2	600 23.62	940 x 850 x 1700 37.01 x 33.46 x 66.93
Ferret Cage	 4P02B700FT	-	2	600 23.62	940 x 845 (+160) x 1700 37.01 x 33.27 (+6.299) x 66.93
Plastic Guinea Pigs Cage	 3W05B100CP	3W05B102CP	5 (Solid base cage type)	280 11.02	860 x 790 x 1700 33.86 x 31.1 x 66.93
	 3W10B100CP	3W10B102CP	10 (Solid base cage type)	280 11.02	1700 x 790 x 1700 66.93 x 31.1 x 66.93
	 3W04B200CP	3W04B202CP	4 (Perforated base cage type)	280 11.02	860 x 790 x 1700 33.86 x 31.1 x 66.93
	 3W08B200CP	3W04B202CP	8 (Perforated base cage type)	280 11.02	1700 x 790 x 1700 66.93 x 31.1 x 66.93

Options

- Ferret Hammock.
- Raised shelf for Guinea Pig Cage.

Accessories

- Nesting box for rabbits.
- Hay feeder for rabbits.

Tecniplast Large Animal Cages are made of EN130S Noryl (GE trademark), guaranteed for 2 years from the date stamp. This plastic is highly resistant, it can be autoclaved at 121°C/250°F for 20 minutes.

Related Products

To maximise operational efficiency and learn more about these products contact your local representative.



900 Series Rack Washer.
Designed specifically to meet the throughput needs of medium to large size facilities, these machines will live up to your highest expectations in terms of maximum cleaning efficiency, safety, environmental impact and quality construction.



650 Series Cabinet Washer.
High throughputs rates, short cycles and environmental friendly.



DS100 Bedding Disposal Station.
Designed to manage and control the levels of contamination generated during the dumping of big rodent cage bedding and animal waste.



"These racks are well designed, very easy to clean and the rabbits are evidently comfortable in them".

Housing
Large Animal Housing
(Rabbits, Guinea Pigs, Ferrets)

- Cost Saving
- Ergonomics
- Flexibility
- Animal Welfare
- Visibility

TECNIPLAST

HEADQUARTERS

ITALY | Tel. +39 0332 809711 • www.tecniplast.it • E-mail: tecnicom@tecniplast.it

AUSTRALIA/NEW ZEALAND | Tel. + 61 2 8845 6500 • www.tecniplast.it • E-mail: info@tecniplast.com.au

FRANCE | Tel. 04 72 52 94 41 • www.tecniplast.fr • E-mail: info@tecniplast.fr

GERMANY | Tel. 08805 921320 • www.tecniplast.de • E-mail: info@tecniplast.de

UNITED KINGDOM | Tel. 0845 0504556 • www.tecniplast.it • E-mail: info@tecniplastuk.com

JAPAN | Tel. + 03 5770 5375 • www.tecniplast.it • E-mail: info@tecniplastjapan.co.jp

USA/CANADA | Toll Free: 877.669.2243 • www.tecniplastusa.com • E-mail: info@tecniplastusa.com



1/ Cost Saving

The components of these cages are easily removable, washable, stackable and autoclavable for good practicality and safety. The Guinea Pig S/S feeder features a plate to prevent the animals from entering and to reduce food spillage.

H Housing Large Animal Housing (Rabbits, Guinea Pigs, Ferrets)

Cost Saving 1



2/ Ergonomics

All the system features external water bottle and feeder. Bottles can be changed and food refilled without disturbing or stressing the animals.

You can extract rabbit and ferret cages without touching the raised shelf, since it features an hooking system with auto-rotation. The rabbit and ferret cages feature a unique quick safe opening/closing system for easy access to the animals.

The Guinea Pig cage feature a S/S front gate door allows the operator to easy access the animals inside the cage without needing to remove the cage itself.

3/ Flexibility

The rabbit and ferret cages feature an interlinking system, allowing to link cages: after removing the side panels, the animals have at their disposal a bigger floor area and can be housed in larger groups for social interaction.



Ergonomics 2



Flexibility 3



Raised shelf

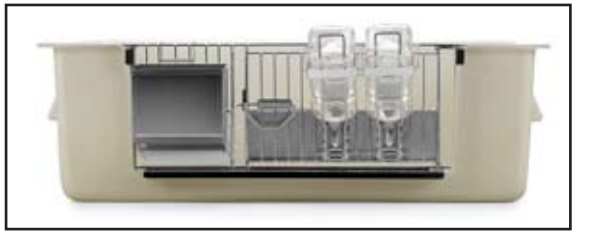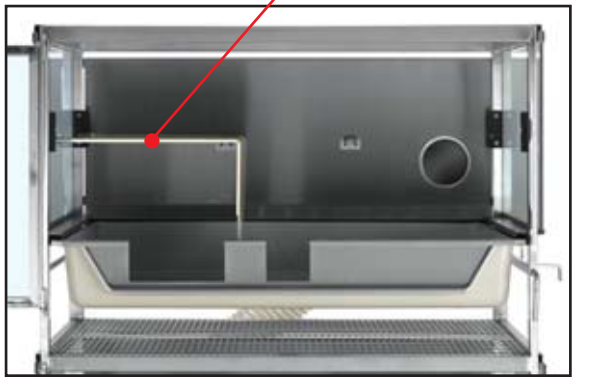

4/ Animal Welfare

All of Tecniplast products for large animal housing comply with US,UK and European Guidelines (ETS 123 revision of appendix A). Rabbit and Guinea Pig perforated Cages are provided with a smooth, perforated floor with rounded corners that provides a warm, quiet surface and cuts out the risks of sore animal hocks.



The rabbit and ferret cages' interlinking system: when the animal is housed in larger groups for social interaction, and this implies natural behavior pattern for the animals.

The ferret cage feature an extractible PVC tube with rigid framework the two cages allowing the animal to move from one cage to the other.

5/ Visibility

Rabbit and Ferret cage feature standard transparent removable panels for good visibility into the cage.

Animal Welfare 4

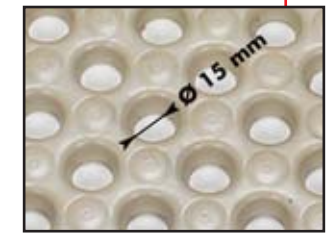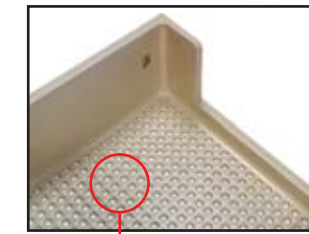

Visibility 5

