


# KEY ACCESSORIES

## LED PLANT GROWTH CHAMBERS

*Models: 7306-22, 7310-22/50/75, 7311-22/50/75, 7312-22/50/75,  
7314-22/50/75, 7315-22, 7317-22/50/75,*

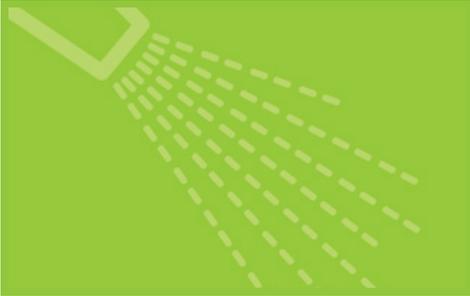

### Controlled Humidity

Precise RH control without messy water pans or flooded unit bottoms. Simple atomizer-based & sensor-controlled system generates no heat and is low maintenance.



### Colored Lighting

4-color LED's for user-adjustable spectra. Includes far-red, red, blue, and white modules. Adjustable by shelf (22 ft<sup>3</sup>) or door (50/75 ft<sup>3</sup>). Available on General and High-Intensity models.



### Condensate Recirculator

Provides a source of Type II equipment water for units with optional controlled RH and eliminates the need for a drain line. Includes sediment & DI filters, and a 7-gallon reservoir. Compact design fits between chambers and rolls out for easy access.



### Chamber Viewing Window

Full-length sliding viewing window makes all plants visible without opening the door and disturbing growth conditions.



### Ramp & Soak Programming

Ramp & soak, 5 linkable profiles, up to 50 steps total. Auto or manual ramp duration, up to 255 process loops. Easy programming, pop-up status screen included.