

EZ-DOCK SMARTBOX® **Euthanasia Station**

Features:



- *Compliant with AVMA Guidelines*
- Cycle starts automatically when cage is docked
- Choose between Adult and Neonates
- LED indicator lights show when cycle is running and when complete
- Automatically shuts off if cage is removed for safety
- Uses house or bottle CO₂ supply at 15-100 psi
- Super E-Z



Scan to watch video



SMARTBOX®
AUTO CO₂ SYSTEM

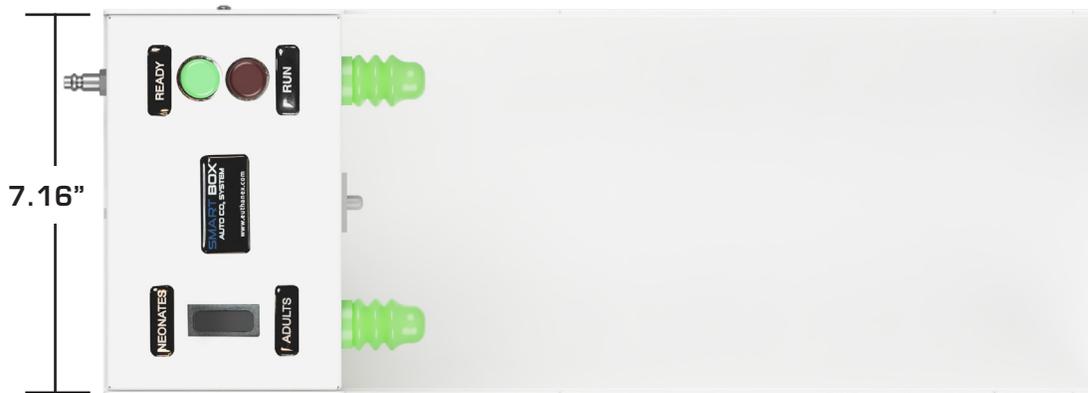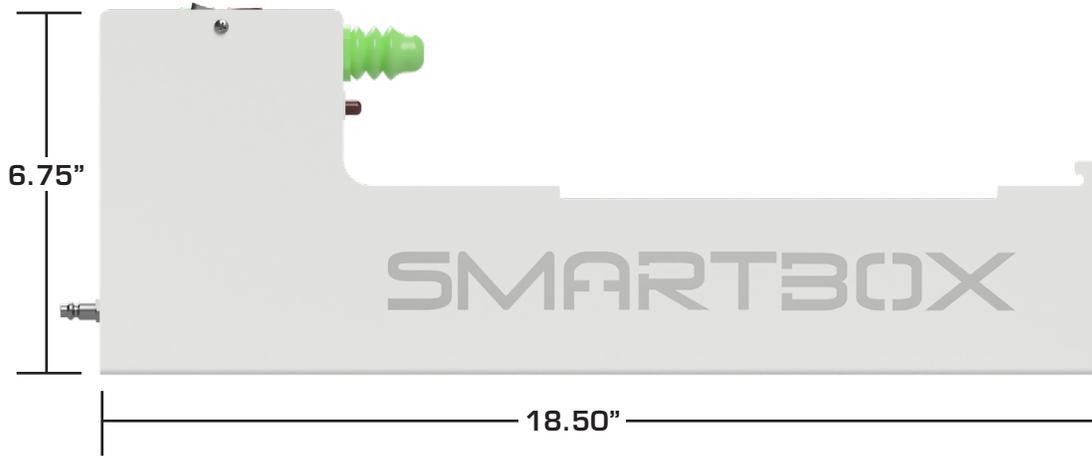

Gas Input Power Input

Power Requirement: 120/240 vac

Incoming Pressure: min 15 psi / max 100 psi

Output flow rate: 3.5L*

Factory Presets: Adult - 5 minute charge

3 minute dwell

Neonate - 5 minute charge

60 minute dwell

*Times and flow rates can be adjusted per request.