



## medical-grade countertop refrigerator

### Features

#### Superior temperature performance

- heavy-duty compressor and forced air cooling
- system maintains +/-1.8 F (1 C) product temperature even with frequent door openings
- quiet compressor and coaxial fans
- non-CFC R134a refrigerant
- no heated defrost cycle required
- dished floor meets NSF requirement

#### Convenience and security

- exterior LED displays product temperature in user-selectable F or C
- "sleep mode" turns off display
- user-programmable, audible and visual high/low alarms
- lockable controller to prevent inadvertent setpoint changes
- remote contacts for connection to central alarm systems
- adjustable refrigeration set points accessed through front controls
- high/low temperature log shows highest and lowest product temperature recorded
- cylinder door lock
- back wall port for third-party probes
- flush door design, compatible with standard automated dispensing system brackets

#### Durability and serviceability

- stainless steel construction on exterior and interior
- rubberized feet
- epoxy-coated wire shelf standard
- heavy-duty edge-mount self-closing hinges
- field-reversible door
- removable back panel for easy service access to compressor, condenser and fans
- dart-type magnetic gasket

#### Warranty

- 2 year parts and labor on refrigeration system
- 5 year parts on compressor

Countertop refrigerators		
Controls	V/Hz/Ph	Item number
LED	115/60/1	REF1

### Options and accessories

Factory-installed keypad locking system (must be specified at time of order)

Wall mount bracket (item # 00969949)

Organizing bins, 6 (item# 00969931)

Extra shelf (item# 00969956)

Temperature surveillance module (item# 00168674) (includes alarm, 6.00" (15.24 cm) chart recorder, remote contacts)

Pedestal, fan-assist, 4.50" (11.43 cm) (item# 01025899)

Pedestal, 12.00" (30.5 cm) (item# 01025881)



Digital data logger kit (item# 01057561)

## Specification

Nominal capacity	1 cu ft (28 L)
W1 Width exterior	18.60" (47.2 cm)
D1 Depth at top	24.00" (61.0 cm)
Depth at bottom	20.00" (50.8 cm)
H1 Height exterior	17.93" (45.5 cm)
Width interior	14.50" (36.8 cm)
Depth interior	10.00" (25.4 cm)
Height interior	11.25" (28.6 cm)
Ventilation clearance with more than 8.00" (20.32 cm) top clearance	2.00" (5.08 cm) for one side
Ventilation clearance with less than 8.00" (20.32 cm) top clearance	6.00" (15.24 cm) for each side
C1 115 V/60/1 electrical	1.8 full load amps, NEMA 5-15P 90° hospital-grade plug low profile, SJT, 18 AWG, 8' (2.4 m) cord
Door	stainless steel
Door lock	cylinder door lock with 2 keys
Door hinges	edge-mount, self-closing
Door swing	standard with self-closing right hinge, field reversible
Storage system	(1) 13.75" x 9.00", (34.9 cm x 22.8 cm) epoxy-coated wire shelf
Gasket	magnetic, dart-style
Temperature display	standard 7 segment LED, default set in F, user-selectable F or C
Display sleep function	button turns off display, press again to turn back on
High and low product temperature alarm	audible and visual, user programmable set points in F or C
Factory set point and operating range	preset at 38 F (3 C) with an adjustable range from 34 F to 55 F (1 C to 13 C)
Alarm	user-programmable high/low; audible and visual alert; audible alarm mute
Alarm delay	user-programmable 0 to 255 min
Automatic defrost	no heated defrost cycle required
Product simulation bottle	(1) 60 ml
Dry contacts	normally open/normally closed (NO/NC) contacts for remote alarming notification
Probes	2 probes : P1 controls refrigeration system; P2 provides product temperature
Refrigerant and charge	R134a refrigerant, 3.4 oz
Energy consumption	1.1 kWh/day
Max heat rejection	190 BTU/hr
Approximate ship weight	95 lb (43 kg)

NOTE: For indoor use only

## Dimensional drawing

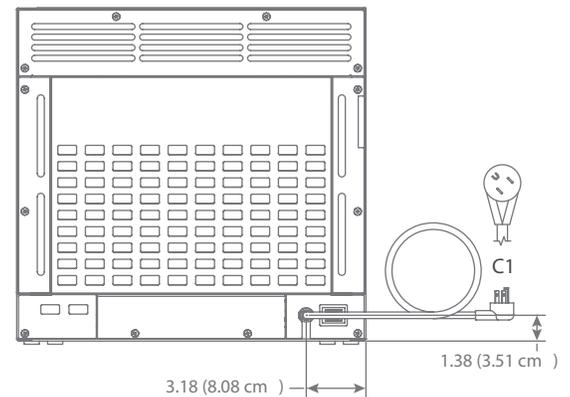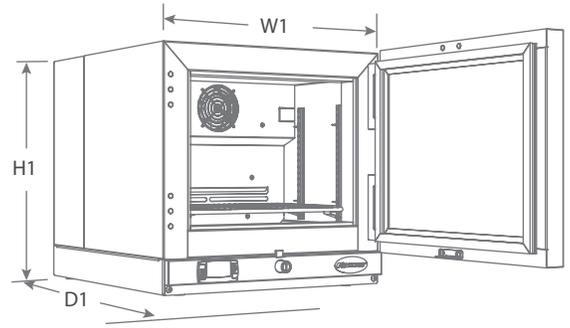

**SHORT FORM SPECIFICATION :** Refrigerator with solid stainless steel door fits on standard 24.00" (61.0 cm) deep countertop and under standard 18.00" (45.7 cm) wall-hung cabinets with 1.0 cubic feet (28 L) of nominal capacity. Includes (1) epoxy-coated wire shelf, adjustable in 1/2" increments. Environmentally responsible R134a forced-air cooled refrigeration system. Bottom-mounted controls display product temperature in user-selectable C or F. Integral high and low temperature alarming. Storage area insulated with CFC-free, high-density polyurethane foam. 8' (2.4 m) power cord with NEMA 5-15P 90° hospital-grade plug. NSF, UL and CUL listed.

Follett is a registered trademark of Follett Corporation, registered in the US. Follett reserves the right to change specifications at any time without obligation. Certifications may vary depending on country of origin.

## REF1 countertop refrigerator

Form #6555 10/14  
© 2014 Follett Corporation



**LABRepCo**  
One Company • Many Solutions

**FOLLETT**