

XC 20 Signature Series

Facio



Features include:


- Engineered specifically for farm use
- Precision manufacturing design

MVE XC 20 Signature	
Maximum Storage Capacity	
Number of Canisters	6
Number of 1/2 cc straws 10/cane	660
Number of 1/2 cc straws 1 Level Bulk	879
Number of 1.2/2.0 mL vials	210
Performance	
LN2 Capacity L	20.5
Static Evaporation Rate L/day*	0.090
Static Holding Time days	227
Normal Working Duration days	142
Unit Dimensions	
Neck Opening <i>in. (mm)</i>	2.16 (55)
Overall Height <i>in. (mm)</i>	26.0 (660)
Outer Diameter <i>in. (mm)</i>	14.5 (368)
Canister Height <i>in. (mm)</i>	11.0 (279)
Canister Diameter <i>in. (mm)</i>	1.65 (41.9)
Weight Empty <i>lb. (kg)</i>	26.0 (11.8)
Weight Full <i>lb. (kg)</i>	62.4 (28.3)

* Static evaporation rate and static holding time are nominal. Actual rate and holding time will be affected by the nature of the container use, atmospheric conditions, and manufacturing tolerances.

TWO Year Standard Warranty
FIVE Year Vacuum Warranty

Tank Features

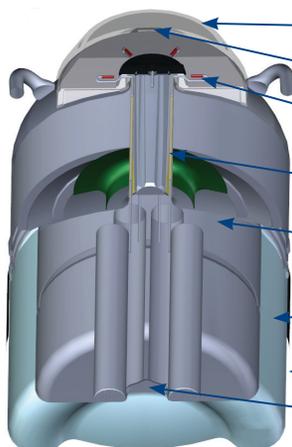

- A** Tamper-proof, hinged lid
- B** Locking tab
- C** Color-coded canister/lid numbering system
- D** High strength neck tube reduces liquid nitrogen loss
- E** Advanced chemical vacuum retention system
- F** Insulation to provide maximum thermal performance
- G** Strong, lightweight aluminum construction
- H** Spider design for easy retrieval and insertion of product canisters

Visit www.chartbiomed.com for more information.