





MVE Doble 22



Features include:


- Engineered specifically for farm use
- Precision manufacturing design
- Tamper-proof locking lid
- Protective shipping container available
- Charges in 2 hours
- Safe for international shipping; IATA Dangerous Good Regulation Exempt (Special Provision A152)

	Doble 22
Maximum Storage Capacity	
Number of Canisters	6 + 1 center
No. of 1/2 cc Straws 10/cane	2400
Vial Capacity	810
Performance	
LN2 Capacity w/o Inventory L est.	20.0
Vapor Capacity L	6.3
Static Evaporation Rate L/day	0.35
Working Duration days for liquid	40
Working Duration days for vapor	18
Unit Dimensions	
Neck Opening in (mm)	3.81 (97)
Overall Height in (mm)	22.0 (559)
Outside Diameter in (mm)	14.5 (368)
Canister Height in (mm)	11 (278)
Canister Diameter in (mm)	3.09 (79)
Weight Empty lb. (kg)	30.0 (13.6)
Weight Charged vapor lb. (kg)	40.0 (18.2)
Weight Full liquid lb. (kg)	68 (30.8)

^{*} Static evaporation rate and static holding time are nominal. Actual rate and holding time will be affected by the nature of the container use, atmospheric conditions, and manufacturing tolerances.

TWO Year Parts Warranty • THREE Year Vacuum WarrantyConforms to MDD 93/42/EEC, the Medical Device Directive for the FU.

Tank Features

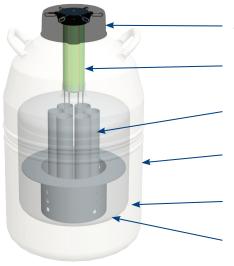

- A Durable composite lid design
- **B** High strength neck tube reduces liquid nitrogen loss
- **C** Canisters
- D Strong, lightweight aluminum construction
- Advanced chemical vacuum retention system
- F Hydrophobic liquid nitrogen absorbent system

Visit www.chartbiomed.com for more information.

Copyright © 2016 Chart Industries Ref ML-A10029 Rev B