


EA-31000 SMARTBOX® Lab Control Unit

SMARTBOX®

This economical unit is designed for small scale use by individual investigators in their laboratories. The Controller automates on/off gas flow at the push of a button, delivering CO₂ to an individual cage or Euthanex Lid. The EA-31000 has a single gas output and includes programmable charge and dwell times. The run and ready lights visually indicate operation and completion.

Operates at a fixed flow rate based on the external CO₂ source.

Compliant with AVMA Guidelines

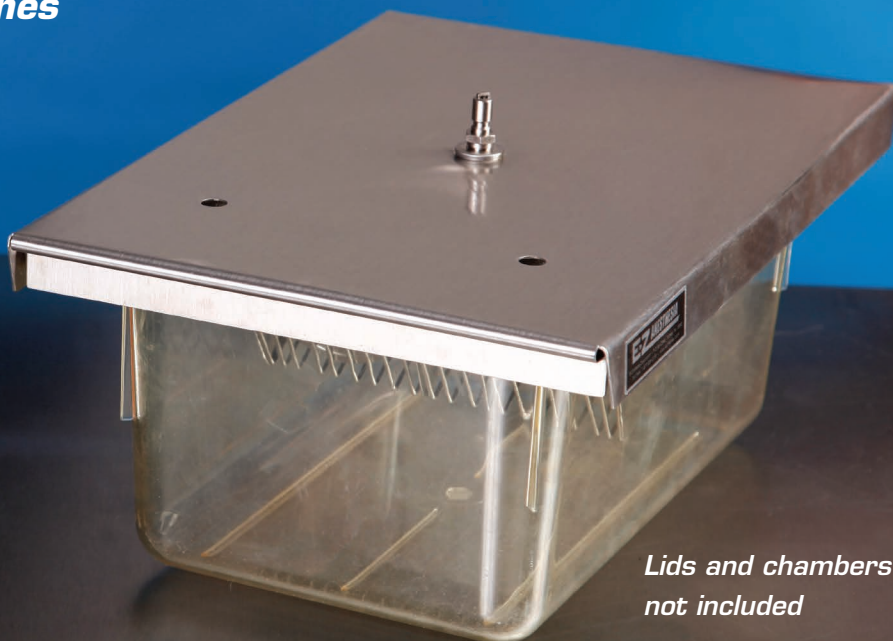


Lids and chambers not included

EA-31000 includes: SMARTBOX® Lab Control Unit and one EA-1130 Lab Unit Hose


Accessories

Lab Control Unit Features:

- Treats animals at individual cage level
- Designed for investigators that need to retrieve animals quickly
- Utilizing the full range of Euthanex Lids or existing induction chamber
- Minimizes operator error and excessive gas use
- Automatic shut-off saves gas and protects personnel
- Lights indicate cycle completion
- Easily attaches to wall

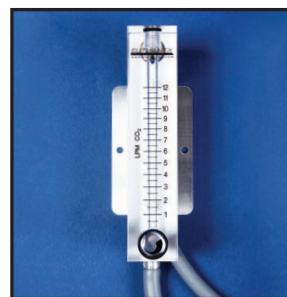

512-A1-BR

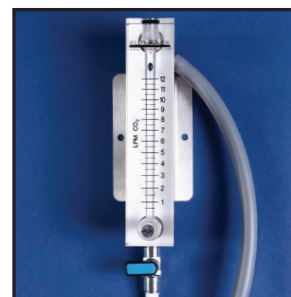

512-T1-BR




EP-PS50




EP-1305


SMARTBOX®

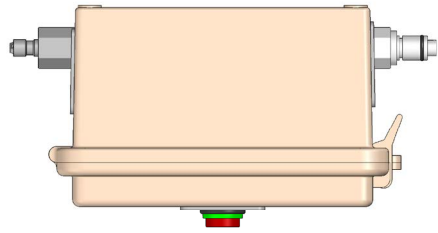


TOP VIEW

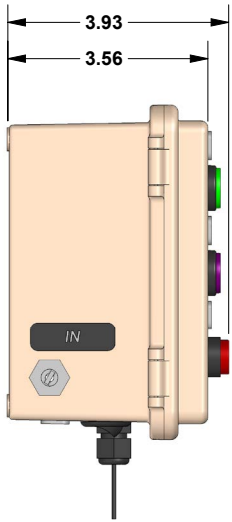

LEFT SIDE VIEW

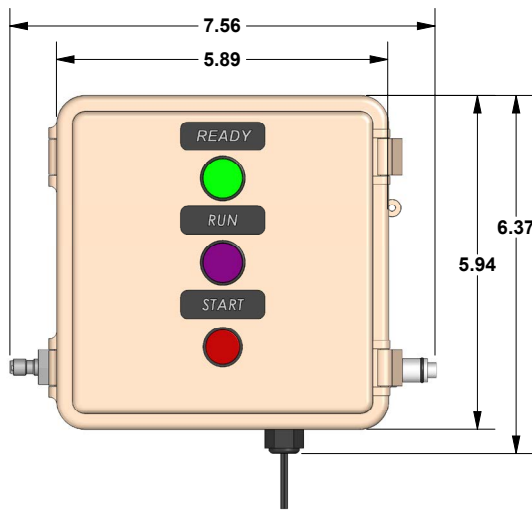


FRONT VIEW

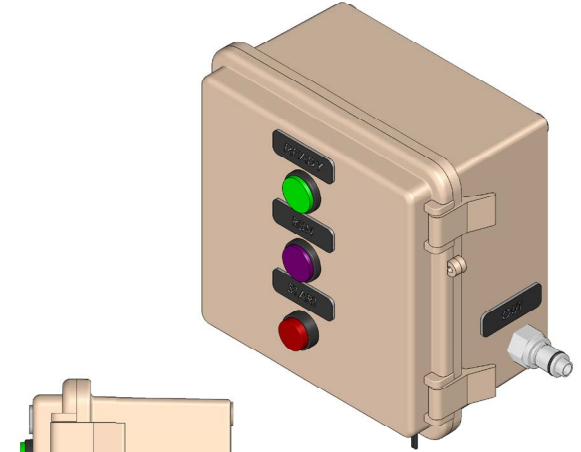

ISOMETRIC VIEW

NOTE:
REQUIRES 110 VAC.

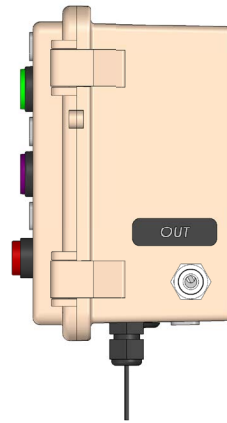

RIGHT SIDE VIEW



BOTTOM VIEW

SMART BOX LAB CONTROL UNIT
EA-31000

E-Z SYSTEMS
www.euthanex.com
1-877-559-0159